

IDEXX **xChek***: Creating a User-Defined Assay for the Pseudorabies Virus gB Antibody ELISA (Filter-Paper Protocol Percent Blocking)

Creating a user-defined assay is an easy process of two tasks—editing an xChek file to allow you to create an assay, and then creating the assay within xChek.

To edit the xChek.ini file:

1. From the task bar, choose **Start > Run** and type “xchek.ini” into the text box. Click the **OK** button.
2. After the xChek.ini file opens, locate the line labeled “AllowNewAssays”.
3. If the value for that line is “0”, change it to “1”. If the value is already “1”, do not do anything.
4. From the menu bar, choose **File > Save**, and then choose **File > Exit**.

To create the assay:

1. From the task bar, choose **Start > Programs > xChek > xChek** to open the xChek program.
2. From the menu bar, choose **Database > Assays > All** to open the Assay dialog box.
3. Click the **New** button and enter “Pseudorabies Virus gB filter paper percent blocking” into the Name field, “PRVgB-fil%” into the Code field, and “Unknown” into the Kit Lot field. Today’s date (mm/dd/yyyy) is automatically entered into the Expiration field.

NOTE: When you begin testing, update the information in the Kit Lot and Expiration fields with the actual kit lot number and expiration date.

The screenshot shows the 'New Assay' dialog box. The fields are filled with the following information:

Name:	Pseudorabies Virus gB filter paper percent blocking
Code:	PRVgB-fil%
Kit Lot:	unknown
Expiration:	6/23/2004
Comment:	

Buttons on the right: OK, Cancel, New, Delete, Print.

Buttons at the bottom: Basic (selected), Calculations, Titers, Bins.

Navigation bar at the bottom: Assays

4. Click the **Basic** button and select or enter the following information for the respective fields:

Case Type: Animal	Samp Filter: 650
Template: SIV Horizontal	Ref. Filter: 0
Species: Cattle	Dilution: 2
Technology: ELISA	Wells: 1

NOTE: Due to a problem with selecting swine as the species type in versions 3.2 and 3.3 of xChek, the species used for this user-defined assay is cattle.

New Assay

Case Type: Samp Filter:

Template: Ref Filter:

Species: Dilution:

Technology: Wells:

5. Click **OK** to save the settings.
6. Click the **Calculations** button and select or enter the following information:

Formulas:

“S” Part of Ratio: None

“P” Part of Ratio: None

“N” Part of Ratio: None

Blocking Factor: $100 * (\text{Neg-Samp}) / \text{Neg}$

First Calculation:

Variable: Blocking

Positive Cutoff: 30

Suspect Cutoff: 30

Comparison: $>=$

New Assay

Formulas

“S” Part of Ratio:

“P” Part of Ratio:

“N” Part of Ratio:

Blocking Factor:

First Calculation

Variable:

Positive Cutoff:

Suspect Cutoff:

Comparison:

Second Calculation

Variable:

Positive Cutoff:

Suspect Cutoff:

Comparison:


7. Click **OK** to save these settings, and then click the **OK** button to save the assay to the database. You can now test for PRVgB (filter-paper protocol percent blocking) using the xChek Assay Management System.

NOTE: xChek does not evaluate controls or results for validity when you use a user-defined assay. You must evaluate the results from each assay in accordance with good laboratory practices. To evaluate your assay validity, please refer to the “Results” section in the insert provided with this test kit.


Setting Up Additional Display Variables

Only OD values are displayed in the reporting options for this user-defined assay. To set up additional display variables, follow the directions below.

1. From the xChek menu bar, choose **Reports > Analyze Cases** to open the Filter Criteria for Analyze dialog box.
2. Enter the desired search criteria and click **OK**. The Select Cases for Analyze dialog box appears.
3. Select the desired cases and click **OK**. The Analyze Report dialog box appears.

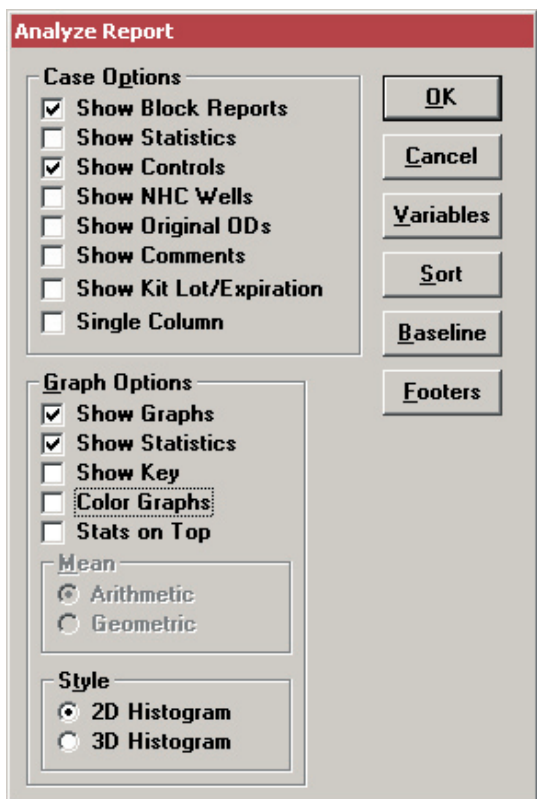

4. Click the **Variables** button to open the Display Variables dialog box.

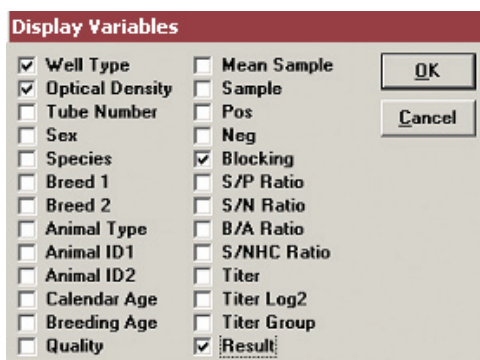

5. Select the desired options and click **OK**. Click **OK** again to save the settings.

For more information, call IDEXX Technical Services at 1-800-548-9997 (in the U.S.) or 1-207-856-0890 (outside of the U.S.), or visit us on the Web at idexx.com.

IDEXX
LABORATORIES

*xChek is a trademark or a registered trademark of IDEXX Laboratories, Inc. in the United States and/or other countries.

© 2004 IDEXX Laboratories, Inc. All rights reserved. • 2794-00

One IDEXX Drive
Westbrook, Maine 04092 USA
idexx.com